

Function

Both the model QVINTETT A and the QVINTETT H give a variety of functions that are easily adjusted with the hand control. Please note that in the sitting position the foot rest goes under horizontal.



Active back in the upright position



Sitting position with foot rest under the horizontal position.



Foot support in the raised position

The beds in the pictures is equipped with bed frame ONYX 4-40.



Central locking wheels with ergonomic pedals.




User-friendly controller.

Active back

The beds back support is active, which means when the back support is raised it lifts upwards and backwards thus reducing the shearing.



Active support system reduces pressure on the patient.

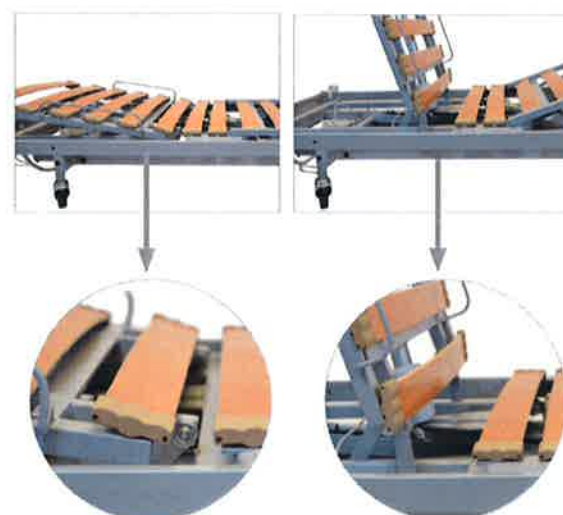

Qvintett H

The foldable bed

QVINTETT H is constructed for easy transport and storage of the bed. The bed frame can easily and simply with out tools is divided into 3 parts, which makes for easy assembly and delivery of the bed. Although it is not necessary to transport the bed divided. Can be delivered just as easy on a transport wagon.



The bed can be transported and stored using a bed transport wagon, a method that works very well for staff.



The bed can be easily divided and packed onto a special bed wagon making it into a unit easily transportable and with a storage area of only 50 x 90 cm with a depth less than 0.5 cm.



Can be transported with a vertical bed lift.



Bed and bed lift in unassembled state.

Weights handling

BED LIFTING	
Lifting support	19 kg
Chassis support	24 kg
Total	43 kg
BED	
Back support	24 kg
Foot support	22 kg
Leg and foot lift	9 kg
Total	55 kg

More product information is available on the back of the brochure.